




J DENTAL CARE

just smile

JD EVOLUTION ^{plus}

Simple insertion, natural stability

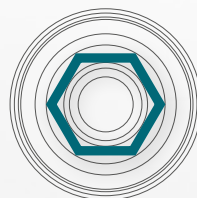

hexagonal connection




Made in Italy

SIMPLE INSERTION, NATURAL STABILITY


features & advantages

- Internal hex connection with lead-in bevel
- Anatomic coronal collar
- Expanding tapered inner body
- Self-cutting capacity in both directions
- Apical cutting blades
- The connection is universally recognized by all the main digital libraries
- Reduces stress concentration on the crestal bone and maximizes bone volume around the implant neck for stable hard and soft tissues
- Makes it possible to achieve high primary stability

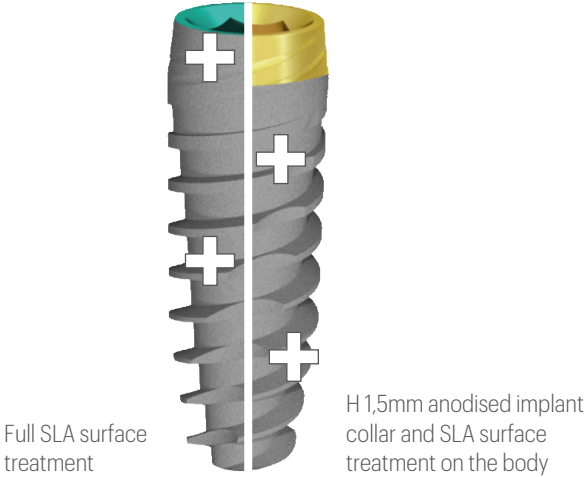

SINGLE PROSTHETIC CONNECTION


Single prosthetic connection for all diameters
Reduced inventory and enhanced working flexibility

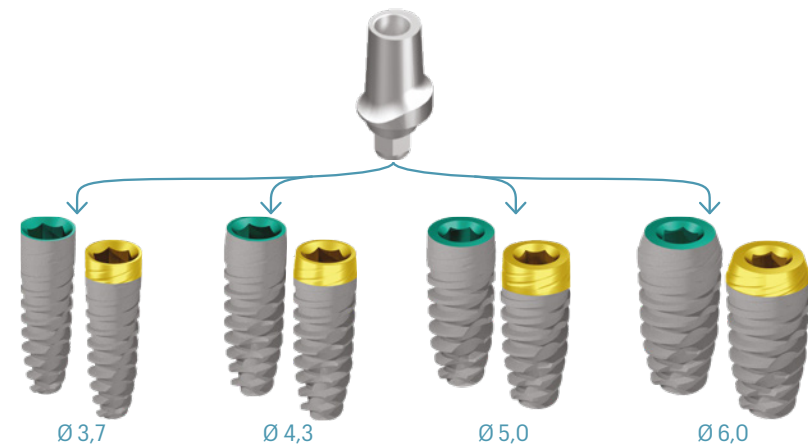

PRODUCT SPECIFICATIONS

The JDEvolution Plus dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	LENGTH	TIP DIAMETER	PLATFORM DIAMETER	ABUTMENT INTERFACE
Ø 3.7	8	2.8	3.7	3.4
	10, 11.5, 13, 15, 18	2.4		
Ø 4.0	8	2.8	3.8	3.4
	10	2.5		
	11.5, 13	2.4		
Ø 4.3	6	3.82	4.0	3.4
	8	3.55		
	10, 11.5, 13, 15, 18	2.8		
Ø 5.0	6	3.85	4.8	3.4
	8, 10, 11.5, 13, 15	3.2		
Ø 6.0	6	4.4	5.0	3.4
	8	4.1		
	10, 11.5, 13, 15	4		

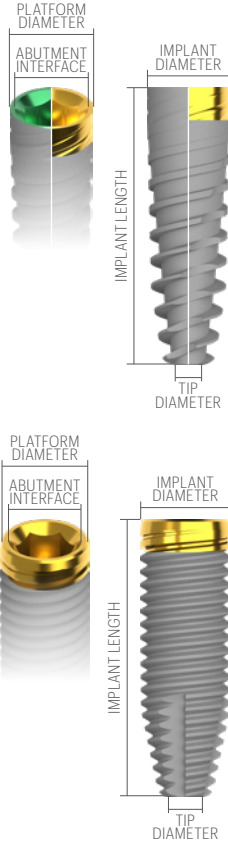

The JDEvolution Plus T dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	LENGTH	TIP DIAMETER	PLATFORM DIAMETER	ABUTMENT INTERFACE
Ø 4.0	6	2.6	3.8	3.4
	8, 10, 11.5, 13, 15, 18	2.4		

Note: All measurements in mm

PRODUCT CATALOGUE

Note: all measurements of diameters and lengths are in mm

Implants:

Ø 3,7

EV37080:	JDEvolution Plus Ø 3.7 L 8
EV37100:	JDEvolution Plus Ø 3.7 L 10
EV37115:	JDEvolution Plus Ø 3.7 L 11.5
EV37130:	JDEvolution Plus Ø 3.7 L 13
EV37150:	JDEvolution Plus Ø 3.7 L 15
EV37180:	JDEvolution Plus Ø 3.7 L 18

Ø 4,0

EV40080:	JDEvolution Plus Ø 4.0 L 8
EV40100:	JDEvolution Plus Ø 4.0 L 10
EV40115:	JDEvolution Plus Ø 4.0 L 11.5
EV40130:	JDEvolution Plus Ø 4.0 L 13

Ø 4,3

EV43060:	JDEvolution Plus Ø 4.3 L 6
EV43080:	JDEvolution Plus Ø 4.3 L 8
EV43100:	JDEvolution Plus Ø 4.3 L 10
EV43115:	JDEvolution Plus Ø 4.3 L 11.5
EV43130:	JDEvolution Plus Ø 4.3 L 13
EV43150:	JDEvolution Plus Ø 4.3 L 15
EV43180:	JDEvolution Plus Ø 4.3 L 18

Ø 5,0

EV50060:	JDEvolution Plus Ø 5.0 L 6
EV50080:	JDEvolution Plus Ø 5.0 L 8
EV50100:	JDEvolution Plus Ø 5.0 L 10
EV50115:	JDEvolution Plus Ø 5.0 L 11.5
EV50130:	JDEvolution Plus Ø 5.0 L 13
EV50150:	JDEvolution Plus Ø 5.0 L 15
EV50180:	JDEvolution Plus Ø 5.0 L 18

Ø 6,0

EV60060:	JDEvolution Plus Ø 6.0 L 6
EV60080:	JDEvolution Plus Ø 6.0 L 8
EV60100:	JDEvolution Plus Ø 6.0 L 10
EV60115:	JDEvolution Plus Ø 6.0 L 11.5
EV60130:	JDEvolution Plus Ø 6.0 L 13
EV60150:	JDEvolution Plus Ø 6.0 L 15



Anodised implant collar:

Ø 3,7

EVA37080:	JDEvolution Plus Anodized Collar Ø 3.7 L 8
EVA37100:	JDEvolution Plus Anodized Collar Ø 3.7 L 10
EVA37115:	JDEvolution Plus Anodized Collar Ø 3.7 L 11.5
EVA37130:	JDEvolution Plus Anodized Collar Ø 3.7 L 13
EVA37150:	JDEvolution Plus Anodized Collar Ø 3.7 L 15
EVA37180:	JDEvolution Plus Anodized Collar Ø 3.7 L 18

Ø 4,3

EVA43060:	JDEvolution Plus Anodized Collar Ø 4.3 L 6
EVA43080:	JDEvolution Plus Anodized Collar Ø 4.3 L 8
EVA43100:	JDEvolution Plus Anodized Collar Ø 4.3 L 10
EVA43115:	JDEvolution Plus Anodized Collar Ø 4.3 L 11.5
EVA43130:	JDEvolution Plus Anodized Collar Ø 4.3 L 13
EVA43150:	JDEvolution Plus Anodized Collar Ø 4.3 L 15
EVA43180:	JDEvolution Plus Anodized Collar Ø 4.3 L 18

Ø 5,0

EVA50060:	JDEvolution Plus Anodized Collar Ø 5.0 L 6
EVA50080:	JDEvolution Plus Anodized Collar Ø 5.0 L 8
EVA50100:	JDEvolution Plus Anodized Collar Ø 5.0 L 10
EVA50115:	JDEvolution Plus Anodized Collar Ø 5.0 L 11.5
EVA50130:	JDEvolution Plus Anodized Collar Ø 5.0 L 13
EVA50150:	JDEvolution Plus Anodized Collar Ø 5.0 L 15
EVA50180:	JDEvolution Plus Anodized Collar Ø 5.0 L 18



Ø 6,0

EVA60060:	JDEvolution Plus Anodized Collar Ø 6.0 L 6
EVA60080:	JDEvolution Plus Anodized Collar Ø 6.0 L 8
EVA60100:	JDEvolution Plus Anodized Collar Ø 6.0 L 10
EVA60115:	JDEvolution Plus Anodized Collar Ø 6.0 L 11.5
EVA60130:	JDEvolution Plus Anodized Collar Ø 6.0 L 13
EVA60150:	JDEvolution Plus Anodized Collar Ø 6.0 L 15



JDEvolution Plus T:

EVT40060:	JDEvolution Plus Ø 4.0 L 6
EVT40080:	JDEvolution Plus Ø 4.0 L 8
EVT40100:	JDEvolution Plus Ø 4.0 L 10
EVT40115:	JDEvolution Plus Ø 4.0 L 11.5
EVT40130:	JDEvolution Plus Ø 4.0 L 13
EVT40150:	JDEvolution Plus Ø 4.0 L 15
EVT40180:	JDEvolution Plus Ø 4.0 L 18



Ø 4,0mm Implants (Tilted):

EVA40150:	JDEvolution Plus Anodized Collar Ø 4.0 L 15
EVA40180:	JDEvolution Plus Anodized Collar Ø 4.0 L 18



Cover screw:

EVCS:	Cover Screw JDEvolution Plus
EVCS6:	Cover Screw Ø 6.0 JDEvolution Plus



Healing Abutment:


EVHA43:	Healing Abutment Ø 4.0 H 3.0 JDEvolution Plus
EVHA45:	Healing Abutment Ø 4.0 H 5.0 JDEvolution Plus
EVHA47:	Healing Abutment Ø 4.0 H 7.0 JDEvolution Plus
EVHA53:	Healing Abutment Ø 5.0 H 3.0 JDEvolution Plus
EVHA55:	Healing Abutment Ø 5.0 H 5.0 JDEvolution Plus
EVHA57:	Healing Abutment Ø 5.0 H 7.0 JDEvolution Plus
EVHA59:	Healing Abutment Ø 5.0 H 9.0 JDEvolution Plus
EVHA63:	Healing Abutment Ø 6.0 H 3.0 JDEvolution Plus
EVHA65:	Healing Abutment Ø 6.0 H 5.0 JDEvolution Plus
EVHA67:	Healing Abutment Ø 6.0 H 7.0 JDEvolution Plus
EVHA753:	Healing Abutment Ø 7.5 H 3.0 JDEvolution Plus
EVHA755:	Healing Abutment Ø 7.5 H 5.0 JDEvolution Plus
EVHA54TS:	Healing Abutment Tissue Ø 5.0 H 4.0 JDEvolution Plus
EVHA56TS:	Healing Abutment Tissue Ø 5.0 H 6.0 JDEvolution Plus



Impression coping:

OPEN TRAY

EVICOT4C:	Impression Coping Open Tray Ø 4.0 JDEvolution Plus
EVICOT5C:	Impression Coping Open Tray Ø 5.0 JDEvolution Plus
EVICOT6C:	Impression Coping Open Tray Ø 6.0 JDEvolution Plus
EVICOTNEC:	Impression Coping Open Tray Non Engaging JDEvolution Plus
EVICOTS:	Impression Coping Open Tray Screw JDEvolution Plus
EVICOT5TSC:	Impression Coping Tissue Open Tray Ø 5.0 JDEvolution Plus

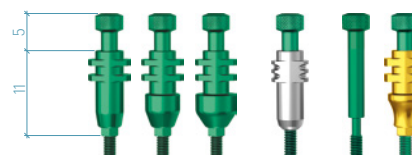

CLOSED TRAY


EVICCT4C:	Impression Coping Closed Tray Ø 4.0 JDEvolution Plus
EVICCT5C:	Impression Coping Closed Tray Ø 5.0 JDEvolution Plus
EVICCT6C:	Impression Coping Closed Tray Ø 6.0 JDEvolution Plus
EVICCTS:	Impression Coping Closed Tray Screw JDEvolution Plus
EVICCT1C:	Impression Coping Closed Tray Direct Press Fit JDEvolution Plus



Temporary Abutment:


EVTAE:	Temporary Abutment Engaging JDEvolution Plus
EVTAEPC:	Temporary Abutment Engaging Peek JDEvolution Plus
EVTANEC:	Temporary Abutment Non Engaging JDEvolution Plus
EVTAEETSC:	Temporary Abutment Tissue Engaging JDEvolution Plus
EVTANETSC:	Temporary Abutment Tissue Non Engaging JDEvolution Plus

Torque recommended 30 Ncm

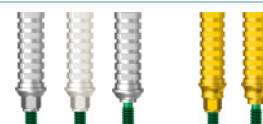


Straight Abutment:

EVNSA4510C:	Straight Abutment Ø 4.5 H 1.5 JDEvolution Plus
EVNSA4520C:	Straight Abutment Ø 4.5 H 2.0 JDEvolution Plus
EVNSA4540C:	Straight Abutment Ø 4.5 H 4.0 JDEvolution Plus
EVNSA5020C:	Straight Abutment Ø 5.0 H 2.0 JDEvolution Plus
EVNSA5040C:	Straight Abutment Ø 5.0 H 4.0 JDEvolution Plus
EVNSA6020C:	Straight Abutment Ø 6.0 H 2.0 JDEvolution Plus
EVNSA6040C:	Straight Abutment Ø 6.0 H 4.0 JDEvolution Plus
EVNSA5020TSC:	Straight Abutment Tissue Ø 5.0 H 2.0 JDEvolution Plus
EVNSA5040TSC:	Straight Abutment Tissue Ø 5.0 H 4.0 JDEvolution Plus

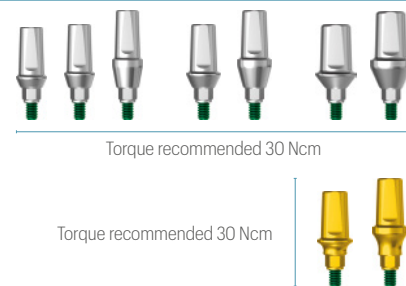


GP Abutment:

EV15GPAEC:	GP Abutment 15° JDEvolution Plus
EV25GPAEC:	GP Abutment 25° JDEvolution Plus
EVGPA40EC:	GP Abutment Engaging Ø 4.0 JDEvolution Plus
EVGPA40NEC:	GP Abutment Non Engaging Ø 4.0 JDEvolution Plus
EVGPAEC:	GP Abutment Engaging Ø 5.0 JDEvolution Plus
EVGPA50NEC:	GP Abutment Non Engaging Ø 5.0 JDEvolution Plus
EVGPAEPC:	GP Abutment Engaging Peek Ø 5.0 JDEvolution Plus
EVGPA60EC:	GP Abutment Engaging Ø 6.0 JDEvolution Plus

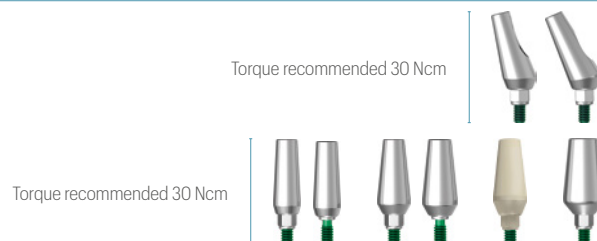


Anatomic Abutment:

EVEA15:	Anatomic Abutment H 1.5 JDEvolution Plus
EVEA30:	Anatomic Abutment H 3.0 JDEvolution Plus
EV15A15:	15° Angulated Anatomic Abutment Ø 5.0 H 1.5 JDEvolution Plus
EV15A15N:	15° Angulated Anatomic Abutment Ø 4.5 H 1.5 JDEvolution Plus
EV15A30:	15° Angulated Anatomic Abutment Ø 5.0 H 3.0 JDEvolution Plus
EV15A30N:	15° Angulated Anatomic Abutment Ø 4.5 H 3.0 JDEvolution Plus
EV25A15:	25° Angulated Anatomic Abutment Ø 5.0 H 1.5 JDEvolution Plus
EV25A30:	25° Angulated Anatomic Abutment Ø 5.0 H 3.0 JDEvolution Plus

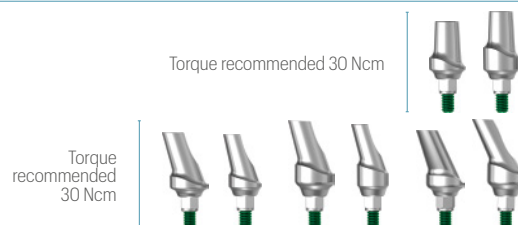


Rapid Abutment:

EVRA15C:	Rapid Abutment H 1.5 JDEvolution Plus
EVRA30C:	Rapid Abutment H 3.0 JDEvolution Plus
EVRA15D:	Rapid Abutment w/o Healing Cap and Coping H 1.5 JDEvolution Plus
EVRA30D:	Rapid Abutment w/o Healing Cap and Coping H 3.0 JDEvolution Plus
EVRAHC:	Rapid Abutment Healing Cap
EVRAIC:	Rapid Abutment Impression Coping
EVRAAR:	Rapid Abutment Replica

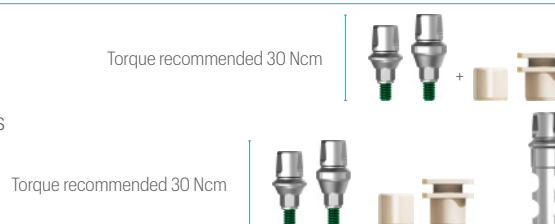


Wax up Abutment:

EVWAEC:	Wax-Up Abutment Engaging JDEvolution Plus
EVWANEC:	Wax-Up Abutment Non Engaging JDEvolution Plus

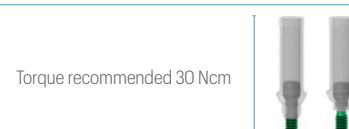


Cobalt-Chrome Abutment:

EVCCEC:	Chrome-Cobalt Abutment Engaging JDEvolution Plus
EVCCNEC:	Chrome-Cobalt Abutment Non Engaging JDEvolution Plus
EVCACCNEC:	Cobalt-Chrome Conical Abutment Non Engaging JDEvolution Plus

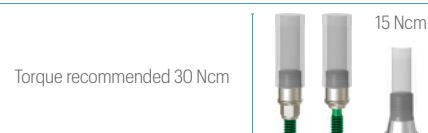

Gold Abutment:

EVGAEC:	Gold Abutment Engaging JDEvolution Plus
EVGANEC:	Gold Abutment Non Engaging JDEvolution Plus

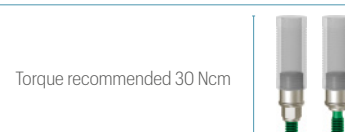

Conical Abutment:

EVNCA10C:	Conical Abutment Straight H 1.0 JDEvolution Plus
EVNCA15C:	Conical Abutment Straight H 1.5 JDEvolution Plus
EVNCA20C:	Conical Abutment Straight H 2.0 JDEvolution Plus
EVNCA30C:	Conical Abutment Straight H 3.0 JDEvolution Plus
EVNCA40C:	Conical Abutment Straight H 4.0 JDEvolution Plus
EVNCA70C:	Conical Abutment Straight H 7.0 JDEvolution Plus
EVNCA90C:	Conical Abutment Straight H 9.0 JDEvolution Plus



EVNCA15TSC:	Conical Abutment Tissue Straight H 1.5 JDEvolution Plus		Torque recommended 30 Ncm
EVNCA30TSC:	Conical Abutment Tissue Straight H 3.0 JDEvolution Plus		Torque recommended 30 Ncm
EVNCA40TSC:	Conical Abutment Tissue Straight H 4.0 JDEvolution Plus		Torque recommended 30 Ncm
EVNCA50TSC:	Conical Abutment Tissue Straight H 5.0 JDEvolution Plus		Torque recommended 30 Ncm
EVCA1725C:	Conical Abutment Angulated 17° H 2.5 JDEvolution Plus		
EVCA1735C:	Conical Abutment Angulated 17° H 3.5 JDEvolution Plus		
EVCA3025C:	Conical Abutment Angulated 30° H 2.5 JDEvolution Plus		
EVCA3035C:	Conical Abutment Angulated 30° H 3.5 JDEvolution Plus		
EVCA1725TSC:	Conical Abutment Tissue Angulated 17° H 2.5 JDEvolution Plus		
EVCA1735TSC:	Conical Abutment Tissue Angulated 17° H 3.5 JDEvolution Plus		
EVA275C:	Conical Abutment Tissue Angulated 17° H 5.0 JDEvolution Plus		
EVCA3025TSC:	Conical Abutment Tissue Angulated 30° H 2.5 JDEvolution Plus		
EVCA3035TSC:	Conical Abutment Tissue Angulated 30° H 3.5 JDEvolution Plus		
EVCA3050TSC:	Conical Abutment Tissue Angulated 30° H 5.0 JDEvolution Plus		
EVA137C:	Conical Abutment Straight Engaging H 1.5 JDEvolution Plus		
EVA138C:	Conical Abutment Straight Engaging H 3.0 JDEvolution Plus		
EVS105:	Screw for Conical Abutment Straight Engaging H 1.5 JDEvolution Plus		
EVS106:	Screw for Conical Abutment Straight Engaging H 3.0 JDEvolution Plus		
EVCASA:	Screw for Angulated Conical Abutment JDEvolution Plus		Torque recommended 30 Ncm
EVCAPS:	Prosthetic Screw Conical Abutment JDEvolution Plus		
EVCAPSL:	Prosthetic Screw Conical Abutment Long JDEvolution Plus		Torque recommended 15 Ncm
EVCAPSA	Prosthetic Screw Angulated for Conical Abutment JDEvolution		
EVCAICOTEC:	Conical Abutment Impression Coping Engaging Open Tray JDEvolution Plus		
EVCAICOTC:	Conical Abutment Impression Coping Open Tray JDEvolution Plus		
EVCAICOTLC:	Conical Abutment Impression Coping Open Tray Long JDEvolution Plus		
EVCAICOT02:	Screw for Conical Abutment Impression Coping Open Tray JDEvolution Plus		
EVCAICOT04:	Screw for Conical Abutment Impression Coping Open Tray Long JDEvolution Plus		
EVCAICCTC	Conical Abutment Impression Coping Closed Tray		
EVCAHC:	Healing Cap for Conical Abutment JDEvolution Plus		
EVCAHCL:	Conical Abutment Healing Cap H 6.0 JDEvolution Plus		
EVCAHC9:	Conical Abutment Healing Cap H 9.0 JDEvolution Plus		
EVCAHCB:	Conical Abutment Healing Cap Peek JDEvolution Plus		
EVCAGPAEC:	GP Abutment Engaging for Conical Abutment JDEvolution Plus		
EVCAGPANEC:	GP Abutment Non Engaging for Conical Abutment JDEvolution Plus		Torque recommended 15 Ncm
EVCATANEWSC:	Temporary Abutment Non Engaging Conical Abutment Smooth for Welding JDEvolution Plus		
EVCATANEC:	Temporary Abutment Non Engaging Conical Abutment JDEvolution Plus		
EVCATANEWC:	Temporary Abutment Non Engaging Conical Abutment for Welding JDEvolution Plus		
EVCAWANEC:	Wax-Up Abutment Non Engaging for Conical Abutment JDEvolution Plus		
EVCAAR	Conical Abutment Replica		Torque recommended 15 Ncm

Prosthetic screw:

EVS:	Prosthetic Screw JDEvolution Plus		
EVSA:	Prosthetic Screw Angulated JDEvolution Plus		Torque recommended 30 Ncm

Interfaces:

EVINTEC:	Interface Engaging JDEvolution Plus		15 Ncm
EVINTNEC:	Interface Non Engaging JDEvolution Plus		
70287C:	Interface Non Engaging for Conical Abutment JDEvolution Plus		
			Torque recommended 30 Ncm
			
			

AS Abutment:

EVASCEC:	AS Abutment Engaging JDEvolution Plus		
EVASCNEC:	AS Abutment Non Engaging JDEvolution Plus		Torque recommended 30 Ncm

TiBase, JDSanBody and Implant Replica CAD CAM:

EVSBCCEC:	ScanBody Engaging JDEvolution Plus		
EVCASBCEC:	ScanBody Conical Abutment JDEvolution Plus		
EVCASBCEEC:	ScanBody Engaging Conical Abutment JDEvolution Plus		
EVCASBCESC:	On Top ScanBody Conical Abutment JDEvolution Plus		
JDSGCIOS	ScanBody Full Arch Tracer		

EVHSB53:	On Top ScanBody Engaging Ø 5.0 H 3.0 JDEvolution Plus		
EVHSB55:	On Top ScanBody Engaging Ø 5.0 H 5.0 JDEvolution Plus		
EVHSB57:	On Top ScanBody Engaging Ø 5.0 H 7.0 JDEvolution Plus		
EVSCPC:	ScanPost JDEvolution Plus		
EVTIBC:	TiBase Engaging Ø4.25 H 0.5 C 4.7 JDEvolution Plus	Torque recommended 30 Ncm	
EVTIBNEC:	TiBase Non Engaging Ø4.25 H 0.5 C 4.7 JDEvolution Plus		
EVTIB15C:	TiBase Engaging Ø4.25 H 1.5 C 4.7 JDEvolution Plus		
EVTIB15NEC:	TiBase Non Engaging Ø4.25 H 1.5 C 4.7 JDEvolution Plus		
EVTIB30C:	TiBase Engaging Ø4.25 H 3.0 C 4.7 JDEvolution Plus		
EVTIB30NEC:	TiBase Non Engaging Ø4.25 H 3.0 C 4.7 JDEvolution Plus	Torque recommended 30 Ncm	
EVTIBNC:	TiBase Long Engaging Ø4.4 H 0.5 C 7.7 JDEvolution Plus		
EVTIBNENC:	TiBase Long Non Engaging Ø4.4 H 0.5 C 7.7 JDEvolution Plus		
EVTIB15NC:	TiBase Long Engaging Ø4.4 H 1.5 C 7.7 JDEvolution Plus		
EVTIB15NENC:	TiBase Long Non Engaging Ø4.4 H 1.5 C 7.7 JDEvolution Plus		
EVTIB30NC:	TiBase Long Engaging Ø4.4 H 3.0 C 7.7 JDEvolution Plus	Torque recommended 30 Ncm	
EVTIB30NENC:	TiBase Long Non Engaging Ø4.4 H 3.0 C 7.7 JDEvolution Plus		
EVA150C:	TiBase Long Engaging Ø4.4 H 0.5 C 7.7 for Angulated Channel JDEvolution Plus		
EVA151C:	TiBase New Non Engaging Ø4.4 H 0.5 C 7.7 for Angulated Channel JDEvolution Plus		
EVA152C:	TiBase Long Engaging Ø4.4 H 1.5 C 7.7 for Angulated Channel JDEvolution Plus		
EVA153C:	TiBase New Non Engaging Ø4.4 H 1.5 C 7.7 for Angulated Channel JDEvolution Plus	Torque recommended 30 Ncm	
EVA154C:	TiBase Long Engaging Ø4.4 H 3.0 C 7.7 for Angulated Channel JDEvolution Plus		
EVA155C:	TiBase New Non Engaging Ø4.4 H 3.0 C 7.7 for Angulated Channel JDEvolution Plus		
EVA133C:	Tibase Interface Engaging H 6.0 for Angulated Channel JDEvolution Plus*		
EVA134C:	Tibase Interface Non Engaging H 6.0 for Angulated Channel JDEvolution Plus*		
EVA135C:	Tibase Interface Engaging H 10.0 for Angulated Channel JDEvolution Plus*	Torque recommended 15 Ncm	
EVA136C:	Tibase Interface Non Engaging H 10.0 for Angulated Channel JDEvolution Plus*		
EVCAITEC:	Interface CAD CAM Engaging for Conical Abutment JDEvolution Plus (H 7.5)		
EVA139C:	TiBase Engaging for Conical Abutment JDEvolution Plus (H 5.0)		
EVCAITC:	Interface CAD CAM Non Engaging for Conical Abutment JDEvolution Plus (H 7.3)	Torque recommended 15 Ncm	
EVCAITBC:	TiBase for Conical Abutment JDEvolution Plus (H 5.0)		
EVS110	Prosthetic Screw Conical Abutment Direct		
EVS111	Prosthetic Screw Conical Abutment Direct Long		
EVS112	Prosthetic Screw for Angulated Channel Conical Abutment Direct		
EVANCN:	Implant Replica CAD CAM New JDEvolution Plus	Torque recommended 30 Ncm	
EVCAARCN	Conical Abutment Replica CAD CAM New		
EVPMNC:	Pre-milled Abutment Blank JDEvolution Plus		

*The EVS: prosthetic screw must be used for straight channels and the EVSA: prosthetic screw must be used for angled channels. Prosthetic screws are not included.
Note: please contact info@jdentalcare.com for the digital library for Exocad, 3Shape and Dental Wings software

Conometric Abutment:

EVCP05C:	Conometric Abutment Straight H 0.5 JDEvolution Plus	Torque recommended 30 Ncm	
EVCP10C:	Conometric Abutment Straight H 1.0 JDEvolution Plus		
EVCP20C:	Conometric Abutment Straight H 2.0 JDEvolution Plus		
EVCP30C:	Conometric Abutment Straight H 3.0 JDEvolution Plus		
EVCP40C:	Conometric Abutment Straight H 4.0 JDEvolution Plus	Torque recommended 30 Ncm	
EV7CP20C:	7.5° Angulated Conometric Abutment H 2.0 JDEvolution Plus		
EV15CP20C:	15° Angulated Conometric Abutment H 2.0 JDEvolution Plus		
EV22CP20C:	22.5° Angulated Conometric Abutment H 2.0 JDEvolution Plus		
EVA300C:	30° Angulated Conometric Abutment H 2.0 JDEvolution Plus	Torque recommended 30 Ncm	
EVA285C:	7.5° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA286C:	15° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA287C:	22.5° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA288C:	30° Angulated Conometric Abutment H 3.0 JDEvolution Plus	Torque recommended 30 Ncm	
EVCPS:	Screw for Conometric Abutment JDEvolution Plus		
EVA130C:	Cap With Pin for Intraoral Welding		
EVA131C:	Cap for Partial or Full-Arch Fixed Prosthesis		
EVA132C:	Anti-Rotation Cap for Single Element		
EVCPIC	Conometric Abutment Impression Coping PEEK		
EVCPAN	Conometric Abutment Replica		

Ball Abutment:

EVBA15:	Ball Abutment H 1.5 JDEvolution Plus
EVBA30:	Ball Abutment H 3.0 JDEvolution Plus
EVBA50:	Ball Abutment H 5.0 JDEvolution Plus
EVBAA	Ball Abutment Replica
EVBAHC	Ball Abutment Elastic Retention Cap w/Container
EVBAN	Ball Abutment Retention Cap



Emi Abutment:

EVEMI15:	Emi Abutment H 1.5 JDEvolution Plus
EVEMI30:	Emi Abutment H 3.0 JDEvolution Plus
EVEMI50:	Emi Abutment H 5.0 JDEvolution Plus
EVEMISCK	Smart cap attachment Emi Abutment
EVEMIHC	Emi Abutment Elastic Retention Cap w/ Container
EVEMIN	Emi Abutment Elastic Retention Cap (1.2 kg)
EVEMIH	Emi Abutment Container
EVEMIIC	Emi Abutment Plastic Impression Coping
EVEMIICS	Emi Abutment Steel Impression Coping
EVEMINT	Emi Abutment White Retention Cap (1.8 kg)
EVEMINY	Emi Abutment Yellow Retention Cap (0.6 kg)
EVEMINP	Emi Abutment Purple Retention Cap (2.5 kg)
EVEMINB	Emi Abutment Retention Cap for Laboratory
EVEMIEIT	Extractor-insertion Tool For Caps
EVEMIA	Emi Abutment Replica



Implant replica:

EVAN:	Implant Replica JDEvolution Plus
-------	----------------------------------



Bone Mill:

JDBMNC:	Bone Mill Ø 5.0 and Guide JDEvolution Plus
JDBM6NC:	Bone Mill Ø 6.0 and Guide JDEvolution Plus
JDBMGNN:	Bone Mill Guide JDEvolution Plus



Short Drills:

JDDICS4	Short Drill Ø 4mm L 6 JDIcon Plus
JDDICS4D	Short Drill Ø 4mm L 6 Dense Bone JDIcon Plus
JDDICS5	Short Drill Ø 5mm L 6 JDIcon Plus



Implant and Prosthetic Drivers:



EVID:	Implant Driver JDEvolution Plus	
EVIDL:	Implant Driver Long JDEvolution Plus	
JDID127	Implant Depth Aligner JDEvolution Plus	
EVSDPF15:	Prosthetic Driver for JDTorque L 15 JDEvolution Plus	
EVSDPF20:	Prosthetic Driver for JDTorque L 20 JDEvolution Plus	
EVSDPF25:	Prosthetic Driver for JDTorque L 25 JDEvolution Plus	
EVSDPF35:	Prosthetic Driver for JDTorque L 35 JDEvolution Plus	
EVSDPF60:	Prosthetic Driver for JDTorque L 60 JDEvolution Plus	
EVSDCAF	Conical Abutment Driver for JDTorque	
EVSDP20:	Prosthetic Driver L 20 JDEvolution Plus	
EVSDP25:	Prosthetic Driver L 25 JDEvolution Plus	
EVSDP30:	Prosthetic Driver L 30 JDEvolution Plus	
JDPD105	Prosthetic Driver JDEvolution Plus for Surgical Driver	
EVSDCA	Conical Abutment Drive	
EVSDPF25A	Angulated Screw Driver for JDTorque L 25	
EVSDPF30A	Angulated Screw Driver for JDTorque L 30	
EVSDPF25AV	Angulated Screw Driver for Conical Abutment Direct JDTorque L 25	
EVSDPF30AV	Angulated Screw Driver for Conical Abutment Direct JDTorque L 30	
EVSUD	Surgical Driver JDEvolution	
JDG	JDGuide	
EVGMC	Implant Mounter for JDEvolution	
EVGMLC	Implant Mounter Long JDEvolution	
JDID102	ZAPA Tool	

Direction indicators:


JDDI	Direction Indicator	
JDDIS	Direction Indicator Short	
JDDI17	Direction Indicator 17°	
JDDI30	Direction Indicator 30°	
JDDI45	Direction Indicator 45°	

DRILLING PROTOCOLS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants

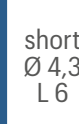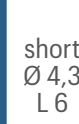	HEALED BONE		POST EXTRACTIVE BONE		  Short Implant L 6mm JDEvolution Plus insertion in maxilla
IMPLANT DIAMETER	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE	
site preparation in maxilla					
Ø 3,7	2,0 2,4 2,8 up to the 1 st laser mark L6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L6mm	2,0 2,4 2,8 at the entrance	2,0 2,4 2,8 at the entrance	
Ø 4,3	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 up to the 1 st laser mark L6-8mm 3,6 at the entrance	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 at the entrance	 short Ø 4,3 L 6 Start the osteotomy with standard twist drill Ø 2.0mm and Ø 2.4mm Complete with the Ø 4mm L 6mm drill code JDDICS4
Ø 5,0	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm	 short Ø 5,0 L 6 Start the osteotomy with standard twist drill Ø 2.0mm and Ø 2.4mm Complete with the Ø 4mm L 6mm drill code JDDICS4D
Ø 6,0	2,0 2,4 2,8 3,2 3,6	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0	2,0 2,4 2,8 3,2 3,6 4,0	
site preparation in mandible					  Short Implant L 6mm JDEvolution Plus insertion in mandible
Ø 3,7	2,0 2,4 2,8 3,2 up to the 2 nd laser mark L 8mm 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 at the entrance	
Ø 4,3	2,0 2,4 2,8 3,2 3,6 up to the 2 nd laser mark L 8mm 4,0 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 up to the 2 nd laser mark L8mm 4,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance	
Ø 5,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 2 nd laser mark L 8mm 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance	 short Ø 4,3 L 6 Start the osteotomy with standard twist drill Ø 2.0mm, Ø 2.4mm and Ø 2.8mm. Complete with the Ø 4mm L 6mm drill code JDDICS4D
Ø 6,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 2 nd laser mark L 8mm 4,8 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4	2,0 2,4 2,8 3,2 3,6 4,0 4,4	 short Ø 5,0 L 6 Start the osteotomy with standard twist drill Ø 2.0mm, Ø 2.4mm and Ø 2.8mm. Complete with the Ø 5mm L 6mm drill code JDDICS5
Ø 6,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 2 nd laser mark L 8mm 4,8 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4	2,0 2,4 2,8 3,2 3,6 4,0 4,4	 short Ø 4,3 L 6 Start the osteotomy with standard twist drill Ø 2.0mm, Ø 2.4mm and Ø 2.8mm. Complete with the Ø 4mm L 6mm drill code JDDICS4D

Note: All measurements in mm

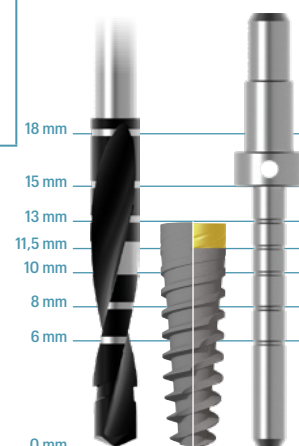

DRILLING PROTOCOLS - Ø 4,0 IMPLANTS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants.

		HEALED BONE		POST EXTRACTIVE BONE	
IMPLANT DIAMETER		SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
site preparation in maxilla					Short Implant L 6mm JDEvolution Plus insertion in maxilla
Ø 4,0	2.0	2.0	2.0	2.0	short Ø 4,0 L 6 Start the osteotomy with standard twist drill Ø 2.0mm and Ø 2.4mm Complete with the Ø 4mm L 6mm drill code JDDICS4
	2.4	2.4	2.4	2.4	
	2.8	2.8	2.8	2.8	
		3.2 up to the 2 nd laser mark L 8mm	2.8 up to the 2 nd laser mark L 8mm		
site preparation in mandible					Short Implant L 6mm JDEvolution Plus insertion in mandible
Ø 4,0	2.0	2.0	2.0	2.0	short Ø 4,0 L 6 Start the osteotomy with standard twist drill Ø 2.0mm, Ø 2.4mm and Ø 2.8mm Complete with the Ø 4mm L 6mm drill code JDDICS4D
	2.4	2.4	2.4	2.4	
	2.8	2.8	2.8	2.8	
	3.2 up to the 2 nd laser mark L 8mm	3.2	3.2 up to the 2 nd laser mark L 8mm	3.2	
	3.6 up to the 2 nd laser mark L 8mm	3.6	3.6 up to the 2 nd laser mark L 8mm	3.6 up to the 2 nd laser mark L 8mm	
		4.0 up to the 1 st laser mark L 6mm			

Note: all measurements in mm



JDENTALCARE

www.jdentalcare.com

JDentalCare, the JDentalCare logo and other trademarks mentioned in this document are, excluding specifically mentioned exceptions, trademarks property of JDentalCare.

The images of the products represented in this brochure are not necessarily to scale.

Note: Not all products manufactured by JDentalCare and represented in this brochure have been regulatory cleared/released for sale in all countries and/or have been licensed in accordance with FDA Rules and Regulations.

For more information about current product assortment and availability, please contact your local JDentalCare distributor. For prescription only.

Caution: Federal (US) law restricts this device to sale by or on the other for a licensed dentist. See instructions for use for more information such as indications, warnings, precaution and contraindications.

REV.58 of 2025-02-07